
# Osaki

## **TP-8500 Massage Chair** **OPERATION INSTRUCTIONS**

**SAVE THESE INSTRUCTIONS**

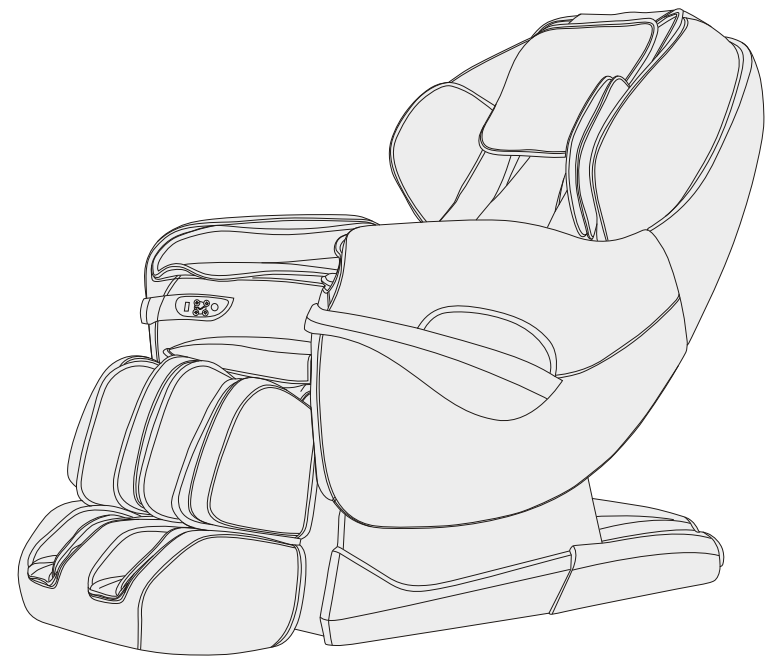Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our Company reserves the right to revise the design and description of this product without any further notice. Please refer to the actual color of this product.

## **CONTENTS**

---

Important Safety Instruction	02
Name and function of components	06
Function	07
Method of usage	08
Product specification	19

## IMPORTANT SAFETY INSTRUCTION

---

**When using an electrical appliance, basic precautions should always be followed, including the following:**

Read all instructions before using (this appliance)

### **DANGER**-To reduce the risk of electric shock:

1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

### **WARNING**-To reduce the risk of burns, fire, electric shock, or injury to Persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use the appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) The appliance is intended for household and indoor use only.
- 14) Connect this appliance to a properly grounded outlet only, see grounding

## Safety And Maintenance

---

instructions .

- 15) Keep children away from extended foot support (or other similar parts)
- 16) Do not use massager in close proximity to loose clothing or jewelry.
- 17) Keep long hair away from massager while in use.

## SAVE THESE INSTRUCTIONS

### **GROUNDING INSTRUCTIONS**

This product must be grounded, if it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock, this product is equipped with a cord having an equipment grounding conductor and a grounding plug, the plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances

### **DANGER**: Improper connection of the equipment grounding

conductor can result in a risk of electric shock, check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded, do not modify the plug provided with the product if it will not fit the outlet, have a proper outlet installed by a qualified electrician

1) This product is for use on nominal 120-volt circuit, and has a grounding plug that looks like the plug illustrated in sketch on the page 13, A temporary adapter that looks like the adapter illustrated in sketches B and C may be used to connect this plug to a 2-pole receptacle as shown in sketch B if a properly grounded outlet is not available the temporary adapter should be used only until a properly grounded outlet (sketch A) can be installed by a qualified electrician, the green colored rigid ear, lug, or the like extending from the adapter must be connected to a permanent ground such as a properly grounded outlet box cover, Whenever the adapter is used, it must be held in place by a metal screw

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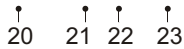
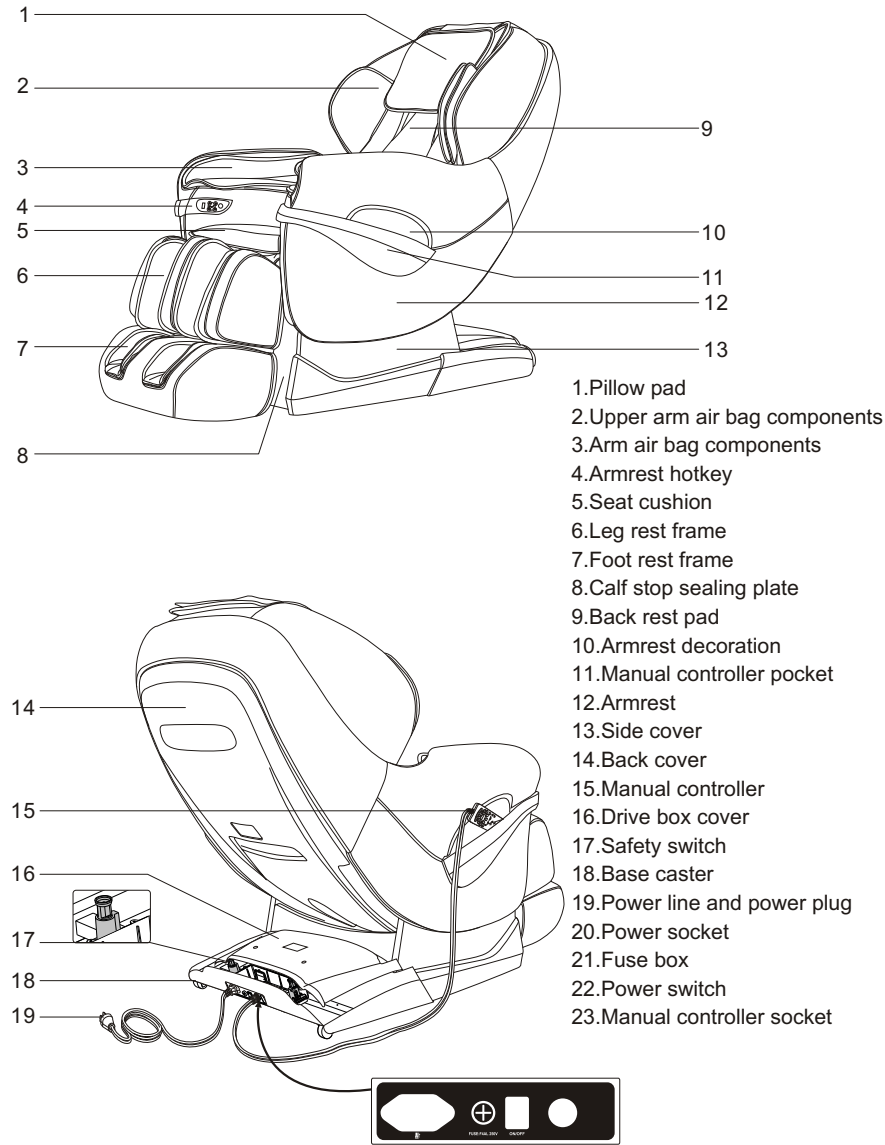
## Function

### Function:

- Designed with a set of vertically movable, four wheel driven muted intelligent massage hands
- Designed with emulation massage methods of shiatsu, kneading, flapping, knocking, simultaneous kneading and flapping
- Set with the shoulder automatic detection and location
- Set with characterized auto massage function: ache-improving, comfort, ease, fast improving experience massage, upper body, lower body six kinds of automatic.
- Set with memory set on function M1 and M2
- Manual massage for the upper body with three options, fixed, partial and overall massage, with five levels of speed adjustable. in the state of flapping, shiatsu and knocking, the width between the two kneading balls is adjustable with three Levels: wide, medium and narrow. On partial and fixed point, the massage hands can move up and down to adjust the massage point accurately.
- Air pressure massage function with three modes: full body air pressure, arms (upper arm set with 4 airbags, arm set with 8 airbags), leg set with 16 airbags, three intensity optional.
- Leg designed with roller massage function at the bottom of foot, two pairs of papillae on toe root rotating massage; on the arch and heel mastoid skin scraping massage.
- Back heating function: carbon fiber far infrared heat source
- Right arm is equipped with a USB charging function and lie-down function.
- The footrest use double stepless extended to suit for different body heights.
- Backrest frame down slides used technology to deviate from the wall by keeping zero space, reach the maximum space saving
- Leg lifting adjusting, the backrest frame lifting adjusting and automatic lay down mode adjustment function
- L shape back track fitting curves of the body, can realize the zero gravity when back down


## Method of usage

Note: The specific displays refer to the actual object of this product!

## Method of usage


### 1. Switch power supply on
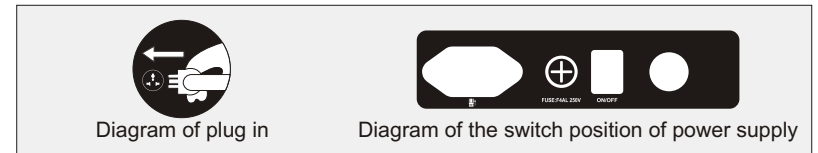
### 2. Start to massage

- 1). Press power button, massage chair automatically reclining, and start massage.
- 2). Automatic shoulder position and massage acupressure point detection, waiting for the completion of testing.
- 3). If there is difference with the massage balls and the actual shoulder position, you can adjust the up and down button to a suitable position. Shoulder position has 10 grades adjustable. It will run automatic massage if there is no adjustment on shoulder position within 10 seconds

Button	Description	Display
Adjust shoulder	Move the massage balls upward to a suitable position	
	Move the massage balls downward to a suitable position	
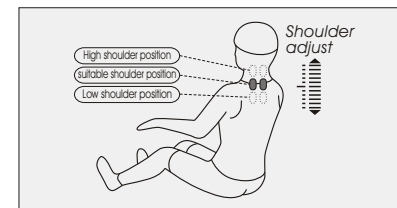
### 4) Choose massage function to massage

Memory function massage

Press button "Upper body / M1" or "Lower body / M2" for 3 seconds, then enter into memory massage function "M1" or "M2", it will revert to the previous memory massage mode.

Function	Button	Display
Memory massage	Upper body / M1	Memory
	Lower body / M2	Memory

## Method of usage

- Characterized auto massage

Function	Button	Feature	Display
Characterized auto massage		Deeply massage on acupressure points, effectively alleviate pain	Auto F1
		Loosen body physique	Auto F2
		Promote blood circulation	Auto F3
		Fast to experience all the massage effect (press for 2S, switch to speed press)	Auto F0
		Massage focus on upper body	Auto F4
		Massage focus on lower body	Auto F5

- Upper body manual massage, lower body manual massage

### 1) Upper body manual massage

Function	Button	Adjust Massage Ball	Display
Technique		Speed And Width	
		Speed	
		Speed And Width	
		Speed	
		Speed And Width	
Speed		Three Grades Speed Adjustable	
Width		Three Grades Widths Adjustable	

## Method of usage

Function	Button	Description	Display	
Back Stretch		Fixed position (fixed point) massage back and forth		
		Small area (partial area) massage back and forth		
		Overall massage back and forth		
		Upper body massage back and forth		
		Lower body massage back and forth		
		Massage point adjust (in fixed point, partial modes you can adjust the massage position)	Massage hand move up	
			Massage hand move down	

- Air pressure function massage

Function	Button	Description	Display
Air pressure mode		Three air pressure modes(full body, arm and leg ) for Lower body	Arm Feet
Intensity		With three grades adjustable	Arm Feet

- Other manual function


Function	Button	Description	Display
Roller		Sole roller: sole roller switch, with 3 speed levels available	
Heating		Far infrared back heating, has warm feeling when it start to work 3 minutes later	

- Other manual function

Function	Button	Description
Calf rise and drop		Calf rest lifting button, press this button, calf restlifting slowly, stop when release it
		Calf rest decline button, press this button, calf rest decline slowly, stop when release it
Backrest rise and drop		Backrest lifting button, press this button, backrest lifting slowly, stop when release it
		Backrest decline button, press this button, backrest decline slowly, stop when release it
Auto lay down mode		with three modes circulation function, when press the button, turn into auto lay downmode.

## Method of usage

### 3.USB recharging function and lie-down function.

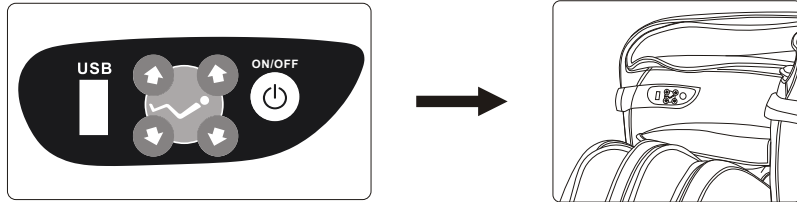
Function	Button	Description
Switch		Press to open and close the massage function.
USB		USB recharging port, recharge smart devices (such as mobile phone).
Calf rise and drop		Calf rest lifting button, press this button, calf rest lifting slowly, stop when release it
		Calf rest decline button, press this button, calf rest decline slowly, stop when release it
Backrest rise and drop		Backrest lifting button, press this button, backrest lifting slowly, stop when release it
		Backrest decline button, press this button, backrest decline slowly, stop when release it

### 4.Turn off power, stop massage

- 1) Press power button during massage progress, all massage function will be stop immediately, the backrest and footrest will return to original position automatically, meanwhile you can set up memory function. After rated working time, all massage functions will be finished, meanwhile you can set up memory function, but backrest and footrest won't return to original position.
- 2) Memory store function: After rated working time or switch off, the hand controller displays "Memory" and glitter, remind the users whether need to store the massage mode before turn off, you can press the memory function buttons "M1" or "M2" as you like, store the massage mode.
- 3) Cut off the power of the whole machine.  
Figure (Switch power off)

## Method of usage
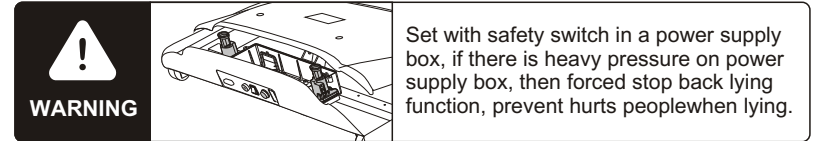
### 5.Armrest disassembly method

Note: this method is only use for cannot easy move into small space, such as too small room door, please do not arbitrarily disassemble.

- Open the master carton, take down the calf rest. (figure 1)
- Remove the back cushion (note: please ensure the heat wire in the cushion was unplugged before remove). (figure 2)

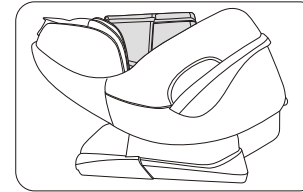
Figure 1

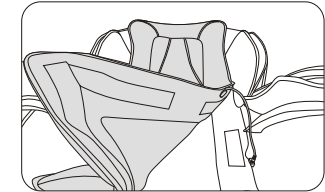
Figure 2

- Zip the backrest kneading cloth open completely, put the cloth on the top of chair (note: please remember to prevent the kneading cloth from dirty or damage). (figure 3)
- The calf rest decoration cover is card slot design, make an effort to pull it when disassembly. (figure 4)

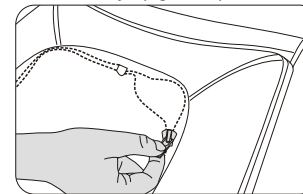
Figure 3

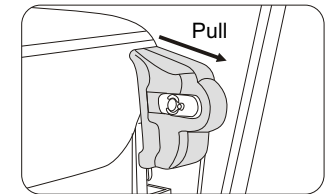
Figure 4

- Unscrew the 8 screws of front dam-board (NEED NOT unscrew the 4 screws on the bottom). (figure 5)
- Open the front dam-board, unplug the joint of air pipe and power cord corresponding to armrest. (figure 6)

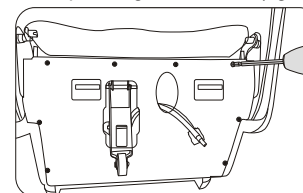
Figure 5

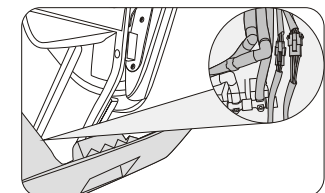
Figure 6

## Method of usage

- Unscrew the screw on the front-end of armrest (note: the screw is in the side of steel frame). (figure 7)
- Unscrew the screw on the front-end of armrest (note: be careful not missing the plastic block for fixed). (figure 8)

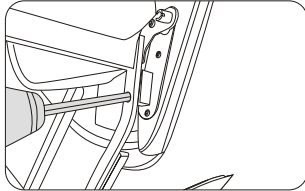
Figure 7

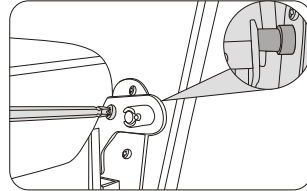
Figure 8

- Unscrew the screw on the back-end of armrest (note: zip the backrest kneading cloth open completely, you may find it). (figure 9)
- Gently lift up the front-end of armrest, pull out the buckle. (figure 10)

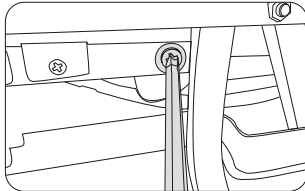
Figure 9

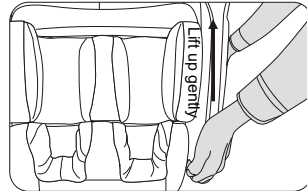
Figure 10

- Gently pull forward the back-end of armrest to pull out the decoration cover, then lift up the back-end of armrest to pull out the buckle. (figure 11)
- Unplug the joint of upper armrest air pipe (figure 12), the armrest removed completely.

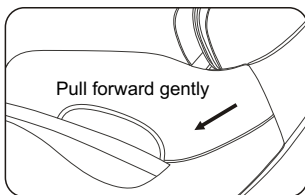
Figure 11

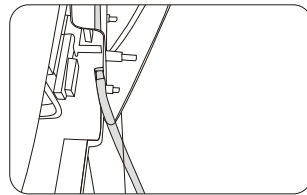
Figure 12

## Method of usage

### 6. Armrest installation method

- Connect the air pipe on the back-end of armrest to the upper arm's. (figure 1)
- Put the air pipe and power cord into the steel frame (avoid the air pipe stuck). (figure 2)

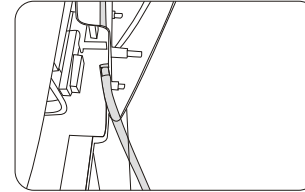
Figure 1

Figure 2

- Hook the back-end of armrest on the steel frame (please note there are two, place it on the upper one). (figure 3)
- Insert the back-end of armrest into the decoration cover of upper arm. (figure 4)

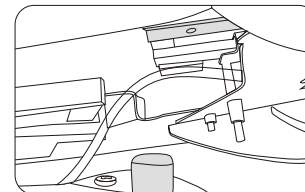
Figure 3

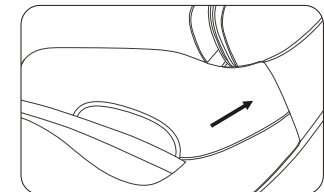
Figure 4

- Hook the front-end of armrest on the steel frame and lock the buckle. (figure 5)
- Gently push the armrest down, double check if the hook well locked. (figure 6)

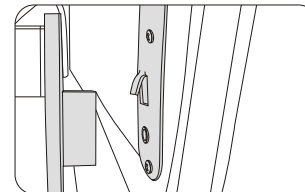
Figure 5

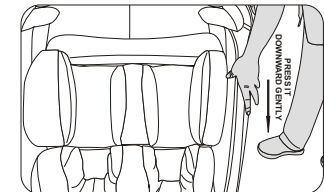
Figure 6

- Put the plastic block between the armrest and calf rest pothook, then tighten the up & down screws. (figure 7)
- Tighten the screws on the back-end of armrest. (figure 8)

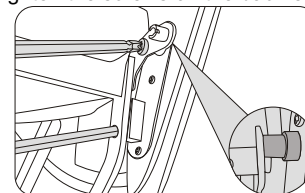
Figure 7

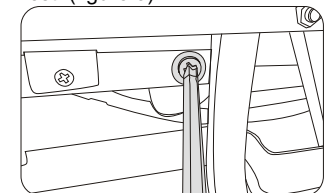
Figure 8

## Method of usage

- Connect the joint of air pipe and power cord (note: there's marks on the air pipe and size of cable connectors is different, you won't miss it). (figure 9)
- Pull the calf rest air pipe and power cord through the hole of front dam-board, then tighten the screws of dam-board. (figure 10)

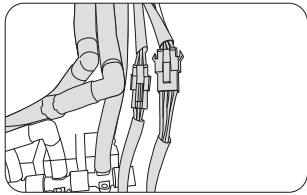
Figure 9

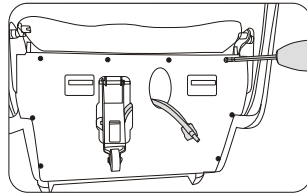
Figure 10

- Zip the backrest kneading cloth up. (figure 11)
- Plug the heat plug of the back cushion into the heat hole on the backrest, put the back cushion tidily on the backrest. (figure 12)

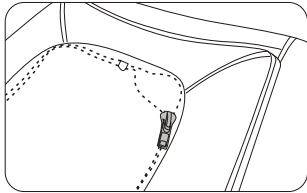
Figure 11

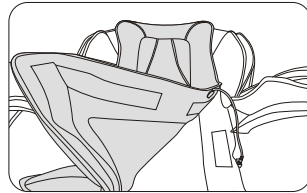
Figure 12

- The calf rest decoration cover is card slot design, make an effort to push it when assembly (note: please keep the cover tight with the side), armrest installation is completed. (figure 13)


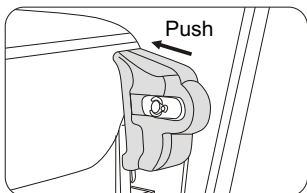
Figure 13


## Method of usage

### 7. Calf installation

- Install the rotating shaft in calf fixed axis first, noticed the hole should be aligned, so that the R shape pin assembly (chart 1)
- After calf is connected with the fixed axis aligned, then insert (chart 2)

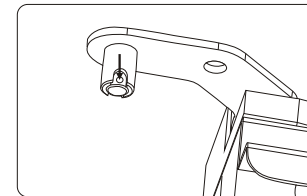
(chart 1)

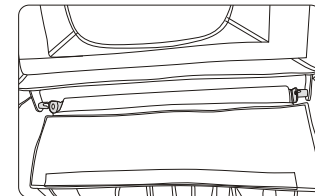
(chart 2)

- Plug R pin, then moved to the other side of the calf, and the plug R pin (chart 3)
- Finished (chart 4)

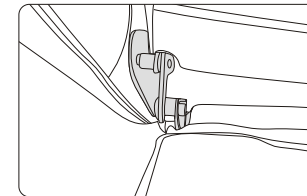
(chart3)

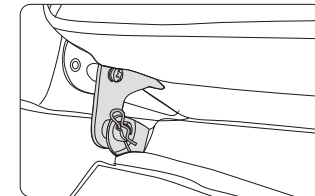
(chart4)

### Note

- When the machine under the working state, with lying and calf rest lifting, you could not apply gravity on the calf, to prevent the massage chair does harm to people before dumping (chart 1).
- Ensure the mechanism hands not in the seat part, otherwise will cause uncomfortable in the Buttock when you seating in the chair (chart 2).

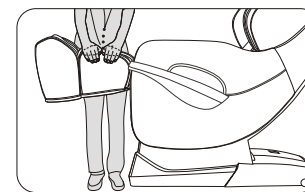(chart1)

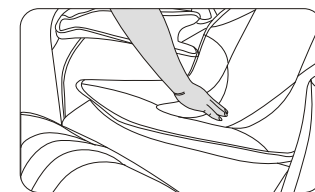
(chart2)

## Method of usage


### 8. Pillow and backrest pad's dismantling and use manual

Using pillows can reduce the kneading massage intensity on shoulder and neck, according to your needs to decide whether the pillow pad is needed (recommend head cushions). Backrest cushion and backrest connected by zipper (1), head cushion and backrest pad connected by thread gluing (2)
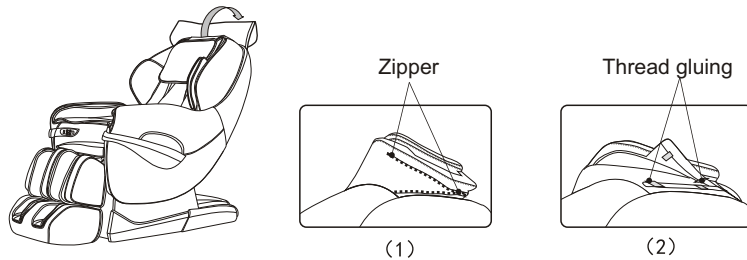
### 9. Floor Protection

Put the heavy massage chairs on the wood floor for a long time may damage the floor, so please place a carpet or other object to avoid this damage.

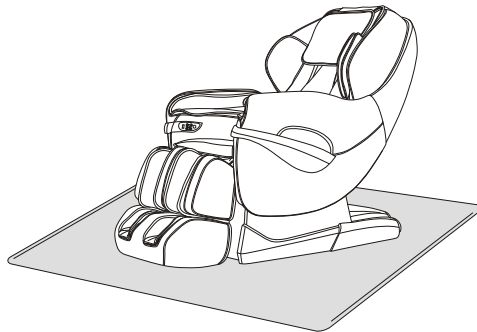
Caution: Do not move the chairs with its wheels on wood floor, rough ground or in an narrow space. Two people are required to lift the chair. During the lifting period, the windscreen and backrest cover grooves as the focal point (Please do not lift armrest directly)


## Method of usage

### 10. Method of movement

Make sure that all wires are far above the ground, tilt the backrest backward to a certain degree (with the gravity center resting on the wheels), push the chair forward or backward with your hands and finally resume the chair to the normal position in a slow and gentle manner.

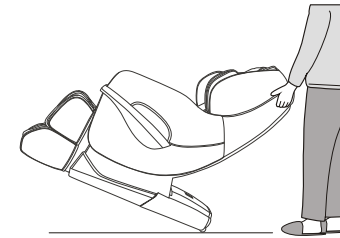Note: During movement you must turn off the power switch first and then unplug the power line and hand controller connecting line
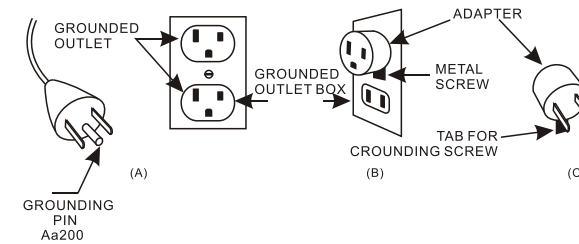
## Product specification

Model: TP-8500

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz

Rated power input: 150W

Rated time: 20Min

Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric & electron parts

MATERIAL	PU	LEATHER	CLOTH	WOOD