

---

KAHUNA MASSAGE CHAIR


# CLOUD







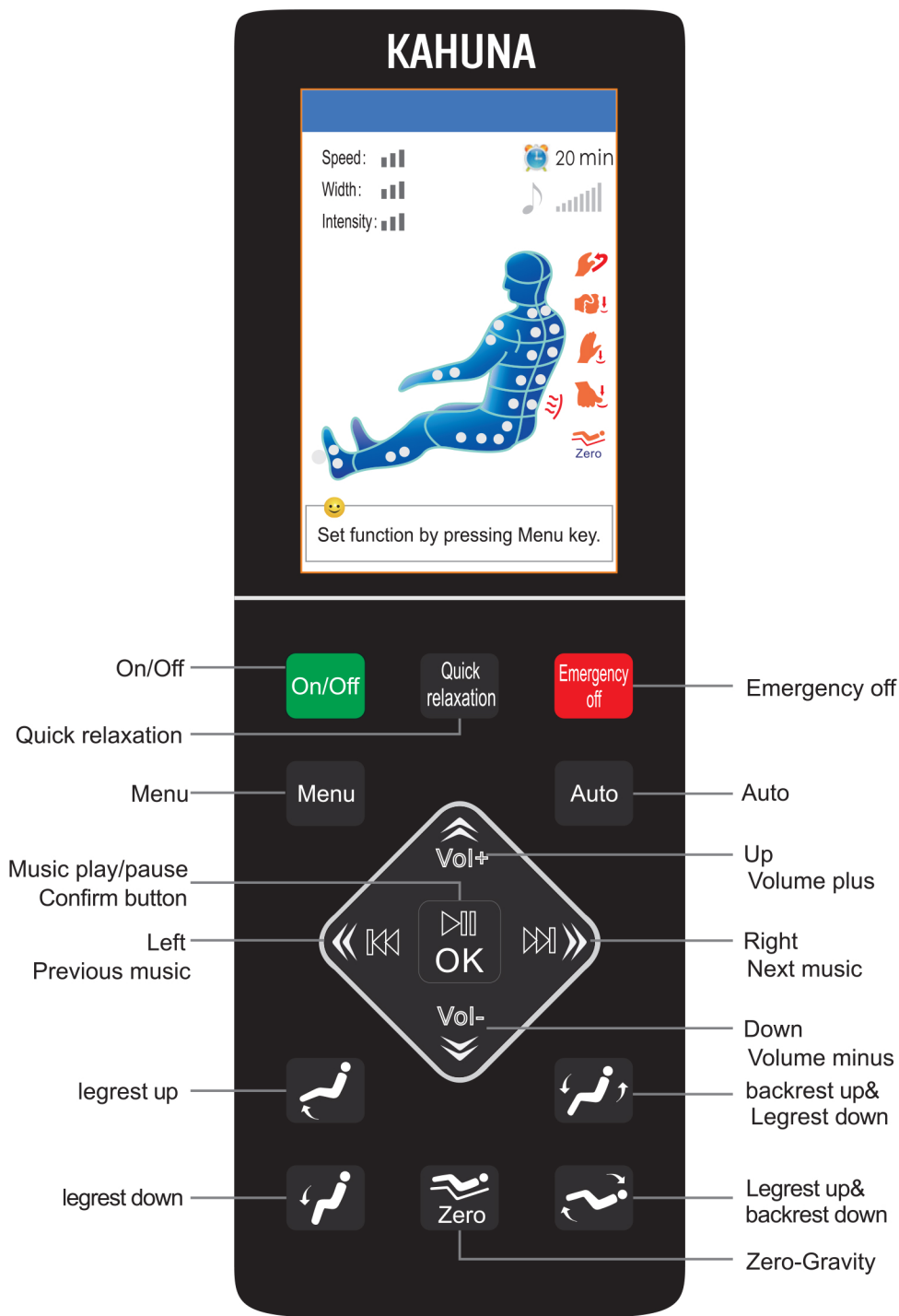


USE AND CARE MANUAL

---




# MANUAL CONTROL




	<p>1. When powered on, the chair will automatically begin in a preparation position. When powered off, the chair will go back to its original position.</p>
	<p>2. In ON status, press it to enter quick relaxation auto program, there is relative status screen display. Default massage time is 8 minutes. Press it again to stop.</p>
	<p>3. Press this button to stop all the on-going massage immediately. All buttons will not work. If user wants to continue massage, please press the On/off button to turn on the controller.</p>
	<p>4. Button in Main interface to switch to menu interface; user can set different functions or programs for massage chair. Press it in Menu again to return to main interface.</p>
	<p>5. This button will help you navigate between several auto and special programs your massage chair has to offer.</p>
	<p>6. Direction button has 4 buttons for up, down, left and right. In Menu interface is the navigator button, user can press up or down to choose menu fast, and press left or right to enter previous or next menu. What's more, up and down button can be used to adjust massage position in acupuncture point tracking detecting finish and point massage. In music operation, left and right button control previous and next song choice; up and down button control volume adjustment. Volume grows one bar every two press.</p>
	<p>7. In music operation, it means to turn on/off the music. In on status and under menu mode, press it to confirm menu choice and enter next menu.</p>
	<p>8. Keep pressing this button, legrest goes down; release it, legrest stop moving. Controller screen will show position adjusting... when adjusting the position of massage chair.</p>
	<p>9. Keep pressing this button, legrest goes up, if you release it, then the legrest stop moving. Controller screen will show position adjusting... when adjusting the position of massage chair.</p>
	<p>10. Keep pressing this button, backrest goes up, legrest goes down, release it, then backrest and legrest stop moving. Controller screen will show position adjusting when adjusting the position of massage chair.</p>
	<p>11. Keep pressing this button, backrest go down, legrest go up, release this button, I backrest and legrest stop moving. When adjusting sitting position, it shows seat position adjusting... in controller screen.</p>
	<p>12. In On status, press this button, massage chair will go into zero-gravity position, press it again to return to default position.</p>

## LCD display graphic description

Main interface

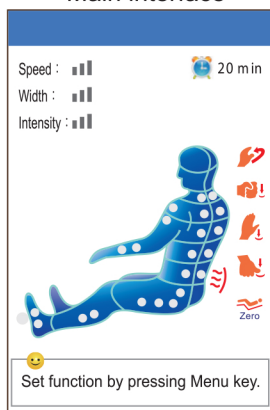

Press "on/off" button when chair power is on, a buzzer beeps once short it. LCD screen display boot screen, after 2-3 seconds it displays the main interface and enter standby mode. Main interface show massage position, method, time, heating current massage status etc., all information will be updated timely when the chair start working.

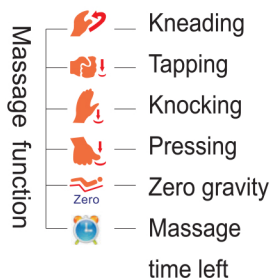

Massage position


- Massage speed (slow, mid, fast)
- Massage width (narrow, mid, wide)
- Air intensity (lowest, low, mid, high, highest)

Menu interface



Press the "menu" button to enter the main interface, it will switch to the menu interface, user can set all massage function. Press the "menu" button in the menu interface will exit and return to the main menu interface.

- Auto massage press the "menu" button, the cursor stop on auto can choose the auto massage mode, and also there are five massage program menu item by default, and 4 special programs.

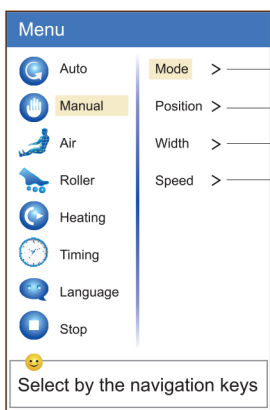

- After Work
- Golfer Mode
- Cloud Squeezes
- Sleep Aid


- Mode >
  - Overall
    - Knead
    - Tap
    - Knock
    - Knead&Tap
    - Shiatsu
  - Narrow
    - Partial
  - Middle
    - Fixed point
  - Low
    - Middle
    - Wide
  - High

Manual press "menu" key, select manual massage by the navigation key, press the "ok" button, there are four options: mode, position, width adjustment, speed adjustment: select the corresponding function, then press the "ok" button to enter the next menu, select the desired massage method, position, width and speed regulating function.


When massage chair in manual fixed point work modes, users can use up/down button to adjust massage position, information for "massage position adjusting..." will be show at the bottom of the display screen.



## LCD display graphic description

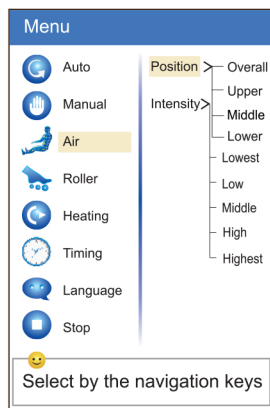


**Air**  
There are two kinds of air pressure related function: Position and Intensity; Choose the corresponding function, then press "ok" to enter next menu and choose required function.

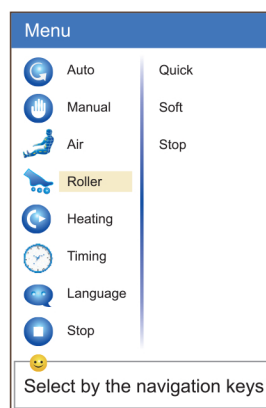

**Roller**  
This menu is for foot massage intensity and off option.



**Heating**  
Health therapy option. The heating symbol will high light display on the interface. With opening and closing two options on the waist.

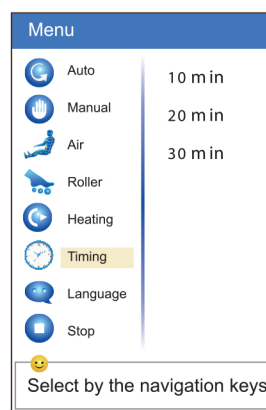

**Timing**  
Time control menu up to 30 min every 10 min interval.



**Language**  
Press "menu" key, select language through the navigation key, then select the preferred language.

**Stop**  
press "menu" key, select stop through the navigation key, it will stop massage, and it will turn back to main interface.




After power is on, the first time to use auto program, massage chair will start acupuncture point position detecting, and "acupuncture point tracking ●●●" will be shown on screen display. After acupuncture point tracking detecting is finished, buzzer will tick 5 times (once every second), meanwhile, the massage machine is

located in shoulders position, if user think the position is not proper, press Up and Down to adjust, "acupuncture point adjusting ●●●" is showing on screen display.

# OPERATION INSTRUCTION AND FUNCTION ILLUSTRATION

## | ON/OFF |

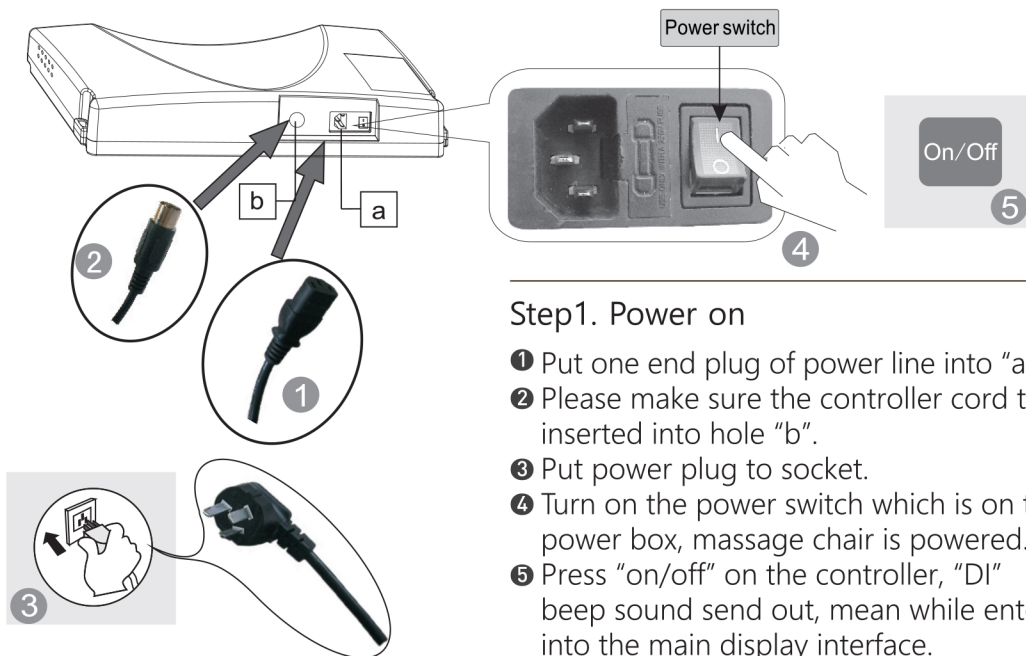


### Step1. Power on

- ① Put one end plug of power line into "a" hole.
- ② Please make sure the controller cord to be inserted into hole "b".
- ③ Put power plug to socket.
- ④ Turn on the power switch which is on the power box, massage chair is powered.
- ⑤ Press "on/off" on the controller, "DI" beep sound send out, mean while enter into the main display interface.



### Step2. Massage

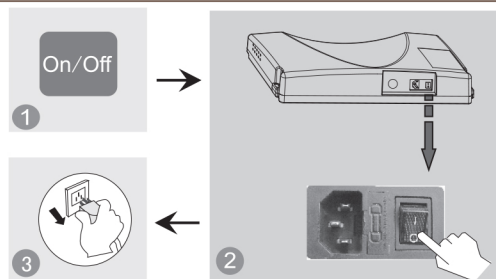- ① Operate according to controller operation instruction.
- ② When massage stops, massage machine will return to the default position, mean while "shutting down" will be displayed on the controller, controller off after 2-3 seconds.



### Step3. Turn off the power

- ① When the chair is on, press on/off button to turn off the massage chair and it will restore automatically.
- ② Turn off the power switch on the bottom of the back of chair after the rollers stay at storage position.
- ③ Unplug the power plug from the socket.

 <p><b>Warning</b></p>	<ol style="list-style-type: none"> <li>① Before using, pull up the pillow, back cushion etc. to check if the PU material is broken, if there is damage, please stop using, unplug the chair, and ask for repairment.</li> <li>② Only plug to proper socket.</li> </ol>	 <p><b>Attention</b></p>	<ol style="list-style-type: none"> <li>① Before sitting down, make sure there is no things pinched between armrest, legrest, backrest and seat of the chair.</li> <li>② Make sure the chair is in returned position before using.</li> <li>③ Do not stand on it when using the massage chair.</li> </ol>
---	--	---	--

# OPERATION INSTRUCTION AND FUNCTION ILLUSTRATION

## Function illustration

### Legrest up

Press this button on controller, legrest goes up. Release the button, legrest stops moving, the angle is locked in this position.



### Legrest down

Press this button on controller, Legrest goes down. Release the button, legrest stops moving, the angle is locked in this position.



### Legrest up and backrest down

Press this button on controller, backrest go down and legrest go up. Release this button, backrest and legrest stop moving, the angle is locked in this position.



### Backrest up and legrest down

Press this button on controller, backrest goes up and legrest go down. Release this button, backrest and legrest stop moving, the angle is locked in this position.



### Zero gravity

Press the button "Zero Gravity", backrest go down and Legrest raise up, the chair go into Zero gravity status automatically. Press the Button again, the backrest and Legrest stop moving, and it will lock at this position. However to change the angels, we can keep comfortable massage status all the time.



### Legrest extend and retract

Legrest can extend or retract within a range of 12 inches.



### Arm massage

Put arms between arm airbags before the embedded airbags before filled with air. It will preform pinch massage techniques to arms and hands after airbags inflating.



### Shoulder massage

Multi-layer extrusion airbags is deployed on both sides of the shoulder. Achieve pressure shoulder stretch massage, relieve shoulder fatigue.



### Foot roller massage

Footrest parts with roller massage and air pressure massage function. Through stimulating the foot reflex zones, and clear the body meridians, regulate and restore human organs function, improve the body resistance and immunity effectively.



- 1 Please make sure there are no people, no pets or any stuff within range for backrest to avoid unexpected hurt or property loss.
- 2 When adjusting the backrest, be careful the gap between the backrest and armrest.
- 3 When the Legrest is moving, please do not stand on or leave from the seat, and do not stop the legrest from moving in case of causing damage to the chair.
- 4 Please refer to control instruction for more details.

## AFTER MESSAGE

### | Restoration of massage chair |

- ❶ Make sure there is no people or pet or anything within the range of backrest and legrest.



- ❷ Make backrest and legrest restore to the normal position by controller.(or press on/off key to make backrest and legrest restore to the normal position)




- ❸ Put the remote control back to the pocket after chair gets reset.




- ❹ Put the pillow pad on to the right position.

- ❺ Switch off the power.

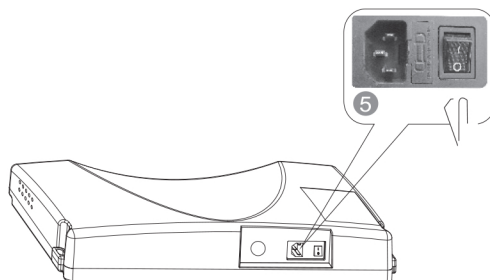

- ❻ Unplug the electricity from the socket.  
(Attention: hold the right plug position)

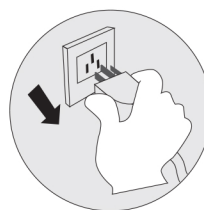

- ❼ Clean the chair if necessary.



### Attention

- After using, press on/off to end all the massage method and wait for restoration.
- Make sure to cut the electricity to avoid any unexpected hurt to children or pets.
- After using, turn off the power, unplug the electricity. Avoid moisture, dirt and aging of insulation in case of causing damage.

## Synthetic leather

### General cleaning

Clean with tidy, soft and dry cloth slightly.



(Chemical or medical cleaning products are forbidden)



Forbidden



### If leather is dirty

- 1 Dip soft cloth into 3-5% neutral detergent solution, squeeze it then clean the stain on the surface gently.
- 2 Wash cloth, squeeze it and wipe with detergent.



- 3 Clean it softly with dry cloth.
- 4 Allow to dry (air blower is prohibited)

## Plastic parts

- 1 Dip cloth into neutral detergent, squeeze it and clean.
- 2 Wipe with water-dipped but squeezed cloth.
- 3 Let it dry fully before using it.



## Clean the controller

- 1 Wipe with soft cloth after squeezing.
- 2 Allow to dry. Propellant, diluent or alcohol are prohibited.



Forbidden



Controller

## Cloth material

- 1 Dip cloth into weak neutral detergent, squeeze it and clean the material surface.

(Diluent, gas, alcohol are prohibited)

Brush the dirty parts with weak neutral detergent. (avoid excessive brushing)

- 2 Excessive cleaning may cause damage to cloth material.
- 3 Wipe with water-dipped but squeezed cloth.
- 4 Allow to dry.



Forbidden



## About cleaning

- Before cleaning, unplug the electricity, don't touch plug with wet hands.
- Make sure hands are dry when plugging to the socket.
- If failed to follow the mentioned instruction, user may get electric shock or burned.





## TECHNICAL DATA

Name	SL-track Massage chair	Power	200 W
Model	CLOUD	NW/GW	Main body : 198 lb/ 246 lb Armrest : 40 lb / 62 lb
Input voltage	AC120V, 60Hz	Packing size	Main body 50 x 29 x 49 inch
Safety design	Class 1		Armrest 44 x 17 x 34 inch



KAHUNA CHAIR PRODUCT WARRANTY REGISTRATION  
<https://www.kahunachair.com/services/>



**Product** Massage Chair

**Model** CLOUD


Developed by AJX Distribution LLC

11125 Knott Ave, Suite D

Cypress, CA 90630

**Tel** 1-855-646-9500

**Email** INFO@KAHUNACHAIR.COM

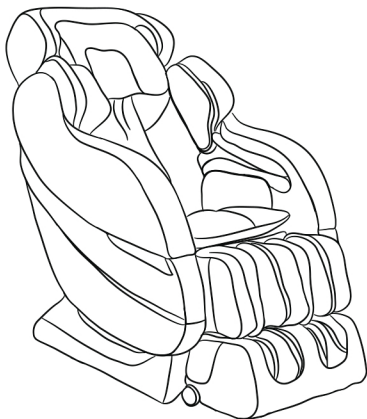

You must register warranty at [www.kahunachair.com/services/](http://www.kahunachair.com/services/)

---

INSTRUCTION MANUAL  
MODEL:CLOUD

Please read this user manual thoroughly prior to operating the massage chair. Store this booklet for future reference.